


FEEDING

TUBE FEEDING

NASOGASTRIC (NG) TUBE

- A NG TUBE IS A FEEDING TUBE THAT IS INSERTED IN ONE SIDE OF THE NARES. THE TUBE GOES DOWN TO YOUR BABY'S STOMACH.
 - SOMETIMES THE FEEDING TUBE MAY BE INSERTED THROUGH THE MOUTH, ESPECIALLY WHEN THEY ARE ON CERTAIN RESPIRATORY SUPPORT, SUCH AS, BUBBLE CPAP
- THE TUBE IS SECURE WITH A SMALL PIECE OF SPECIAL TAPE.
- THE TUBE WILL BE USED ANYTIME THE BABY CANNOT TAKE ALL NECESSARY FEEDINGS BY MOUTH.
- THE FEEDING TUBE CAN BE IN PLACE FOR 7 DAYS AT A TIME.

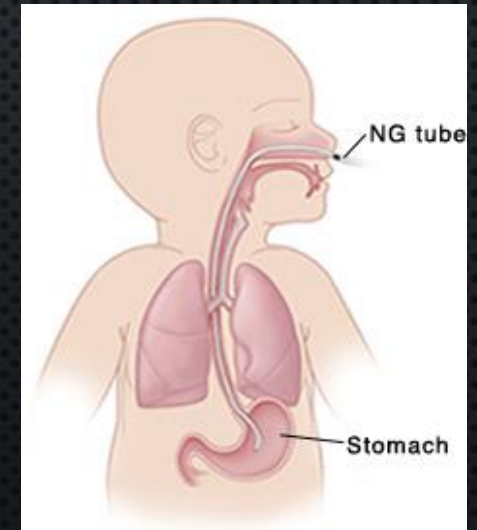

NASOGASTRIC (NG) FEEDING TUBE

- IT IS IMPORTANT TO CHECK THE POSITION OF YOUR BABY'S FEEDING TUBE PRIOR TO USING IT.
 - GOOD PLACEMENT MEANS THE FOOD WILL BE GOING INTO YOUR BABY'S STOMACH AND NOT INTO YOUR BABY'S LUNGS.
 - TO CHECK PLACEMENT OF FEEDING TUBE: LOOK FOR THE NUMBER JUST OUTSIDE OF THE NARES AND CONFIRM WITH YOUR NURSE THAT THIS NUMBER IS APPROPRIATE FOR YOUR BABY'S AGE AND SIZE.

VERIFYING MILK

- IT IS IMPORTANT THAT PRIOR TO EVERY FEEDING THE SYRINGE OR BAG OF MILK IS VERIFIED TO MATCH YOUR BABY'S NAME
- YOU CAN DO THIS BY LOOKING AT THE LABEL ATTACHED TO THE SYRINGE OR BAG
 - NOTIFY YOUR NURSE IF THE LABEL DOES NOT MATCH YOUR BABY'S LAST NAME

HOW TO WARM YOUR BABY'S MILK

- WASH OR GEL HANDS
- USE NEW PENGUIN THERMA-LINER BAG FOR EACH 12 HOUR SHIFT
- OPEN FRONT, BLUE ZIP-LOCK SECTION, AND FILL WITH PLAIN WATER TO THE DOTTED LINE
- ZIP CLOSED
- PLACE THERMA-LINER BAG IN TO THE PENGUIN WARMER
- AFTER VERIFYING NAME ON SYRINGE, PLACE SYRINGE OR BOTTLE OF MILK IN THE TOP ORANGE SECTION OF THE THERMA-LINER BAG



WARMING YOUR BABY'S MILK

- PRESS FAR LEFT “WARM/OFF” BUTTON ONCE
- MACHINE WILL START AUTOMATICALLY FOR <10 ML VOLUMES OF MILK
- PRESS FAR RIGHT “SELECTION” BUTTON ONCE FOR 10-19 ML, TWICE FOR 20-29 ML, THREE TIMES FOR 30-39 ML VOLUMES, ETC.
- RED “WARMING” LIGHT WILL ILLUMINATE AND MACHINE WILL START
- GREEN “READY” LIGHT WILL ILLUMINATE WHEN MILK IS WARMED
- PRESS FAR LEFT “WARM/OFF” BUTTON TO TURN OFF THE WARMER

USING MEDELA SYRINGE PUMP

- [HERE IS A LINK TO A SHORT VIDEO ON THE MEDELA SYRINGE PUMP: RIGHT CLICK ON THE LINK BELOW AND SELECT "OPEN HYPERLINK"](#)
- [HTTPS://WWW.MEDELABREASTFEEDINGUS.COM/VIDEO-PAGE/112/ENTERAL-FEEDING-PUMP-IN-SERVICE](https://www.medelabreastfeedingus.com/video-page/112/enteral-feeding-pump-in-service)



MEDELA SYRINGE PUMP

- WASH OR GEL HANDS
- SQUEEZE BLACK CLUTCH LEVER TOWARD WHITE TAB TO RELEASE CLUTCH
- WHILE SQUEEZING BLACK TAB SLIDE SYRINGE DRIVER ALL THE WAY TO THE RIGHT
- RELEASE CLUTCH LEVER
- LIFT CLEAR SYRINGE CLAMP AT CENTER
- INSERT SYRINGE WITH TUBING ATTACHED AND LOWER CLEAR SYRINGE CLAMP MAKING SURE THE SYRINGE TAB IS LOCKED INTO THE METAL TAB
- SQUEEZE BLACK CLUTCH LEVER AND SLIDE SYRINGE DRIVER TO THE LEFT UNTIL IT CONTACTS THE END OF THE SYRINGE PLUNGER

MEDELA SYRINGE PUMP

- PRESS POWER ON
- PRESS #1 (VOL/TIME)
- PRESS ENTER/YES
- SYRINGE IS NOW LOADED
- PRESS ENTER/YES AGAIN
- SYRINGE SIZE IS DETECTED

ENTERING VOLUME TO BE DELIVERED

- NEXT YOU WILL PUT IN A VOLUME OF MILK:
 - PRESS 2 AND 0 FOR 20 ML OF MILK
 - PRESS 3 AND 0 FOR 30 ML OF MILK ETC.

ENTERING TIME FOR FEEDING

- PRESS ENTER/YES AGAIN
- NEXT YOU WILL PROGRAM A TIME
 - PRESS 3 AND 0 FOR 30 MINUTES
 - PRESS 6 AND 0 FOR 60 MINUTES
 - PRESS 1 AND 30 FOR A 90 MINUTE FEEDING

- PRESS ENTER/YES AGAIN
- PRESS DELIVER
- CONNECT TUBING FROM SYRINGE TO BABY'S FEEDING TUBE
- PRESS ENTER/YES AGAIN
- A GREEN LIGHT WILL FLASH
 - THIS IS HOW YOU KNOW THE FEEDING IS GOING
 - IF THERE IS A RED FLASHING LIGHT, THE FEEDING DID NOT START

IN ----- ON ----- CONNECT ----- GREEN